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
2. Sauger montieren

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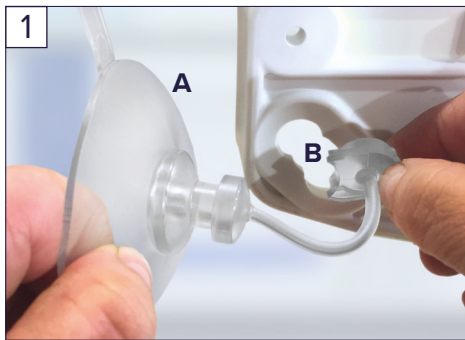
Quick instructions – initial installation

1. Remove from box

1. Lift back rest

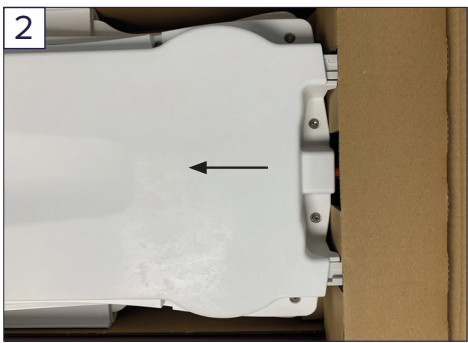
2. ... and remove.

3. Remove seat unit.

2. Mount suction cups

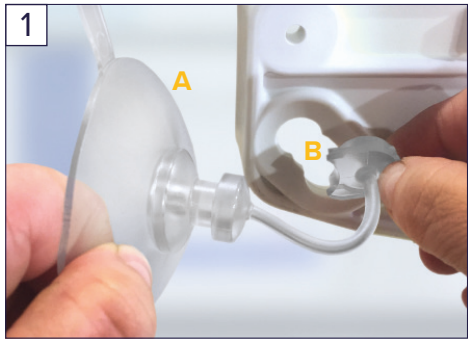
1. To insert the rear suction cups, begin by first inserting the movable plugs (B) into the underside of the base plate.

2. Push the mushroom head (A) of the suction cup toward the small opening by gently rotating it into its final position. Make sure the vent tabs (C) are pointing downward.


3. Push the plug from the upper side into the plate opening until it is flush.

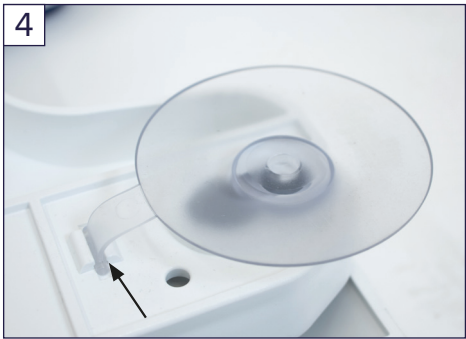4. Insert the front suction cup in the same way. In this case, however, you must also insert the vent tabs into the base plate slots from the side. This will automatically vent the front suction cups during removal.


3. Insert back rest into seat unit

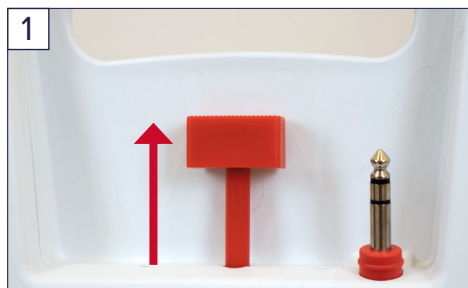
1. Pull the red control lever to the uppermost position and hold it there throughout the entire process of inserting the backrest into the seat surface.


2. To insert the back rest, position the control lever mounting fixture over the bolt in the seat unit.

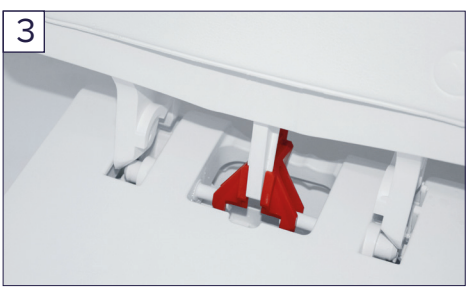3. Insert the drive unit with the integrated back rest into the seat unit as shown in Fig. 2. Now reposition the drive unit down and to the back.

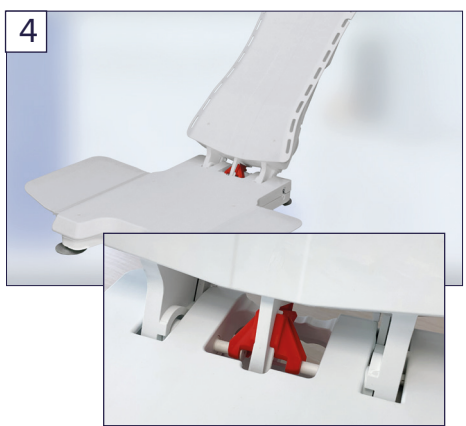
4. Charging the battery • Plugging in the manual control • Care and maintenance

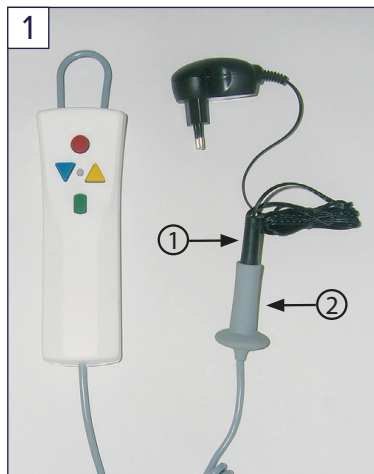Charge the battery completely before using for the first time.

1. First plug the power supply adapter 1) into the socket 2) attached to the manual control and plug the power supply adapter into an electrical outlet. The control light will be red while charging and will turn green when the battery is charged. Once the battery is fully charged, disconnect the manual control from the charger unit. First remove the power plug from the electrical outlet.

Storage at roughly 15–20 degrees Celsius is recommended. Do not place the manual control on a source of heat. Avoid direct sunlight.

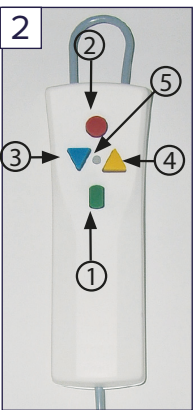2. Manual control description:

- 1. Green: On button
- 2. Red: EMERGENCY OFF button
- 3. Blue: Down button
- 4. Yellow: Up button
- 5. Control light

Incorrect!

Correct!

3. Plug in the manual control behind the back rest. To do this, take the socket attached to the coil cord and slide it onto the plug until it noticeably clicks into place and the red ring on the headphone-style plug is no longer visible.

Do not clean using any highly concentrated chlorine-based cleaning agents!